

HAMILTON

Verso® Dx

Automated Storage Solution for Post-Analytics



Revolutionize Your Post-Analytics

The innovative Verso Dx is an automated post-analytic storage system that caters to different workflows with its capability to load samples either manually or via a lab track. The flexible system combines the cold room and storage retrieval module within one single platform.

Verso Dx can be designed and built in different sizes and features a storage capacity ranging from 30,000 to 1 million tubes. The system can store tubes at ambient or the common temperature of 4°C, while also providing the ability to store samples at -20°C for special applications.

Verso Dx users benefit from streamlined specimen handling thanks to the system's automated retrieval and discard functions. The advantages compared to manual tube storage are amplified when Verso Dx is connected to a laboratory information system (LIS), allowing users to retrieve jobs processed directly through an LIS request. With the system's versatility, usage can range from classical post-analytic storage to other storage and sorting tasks in the laboratory workflow.

Flexible system usage:

- Automated High-Speed Discard
- Manual Input/Output
- Input/Output through TLA



Setting the Green Standard in Cold Storage

Hamilton's full line of automated storage systems is now available with natural refrigerants. The full transition reinforces our commitment to sustainable innovation while helping customers reduce energy consumption, support ESG goals, and future-proof operations without sacrificing sample integrity or system performance.

Maximize Efficiency

Verso Dx allows labs to combine a storage retrieval module and the traditional cold room into one fully automated specimen storage system. The large Input/Output module can be manually loaded with up to 1,300 sample tubes at a time, though a direct connection to TLA tracks is also available.


Automated sorting, retrieval, and discard are integrated, allowing labs to streamline their post-analytical sample management. With a Verso Dx in your lab, you can maximize the efficiency of your post-analytic workflows and give lab personnel more time to focus on high-value tasks.



**LEARN
MORE**

Flexible System Configuration

- Storage capacity of 30,000 to 1 million tubes
- Integration compatibility with major TLA providers
- System can be designed and built to fit your lab's space requirements
- Available storage temperatures:
 - Ambient
 - +4°C
 - -20°C
- Redundant refrigeration for superior sample safety

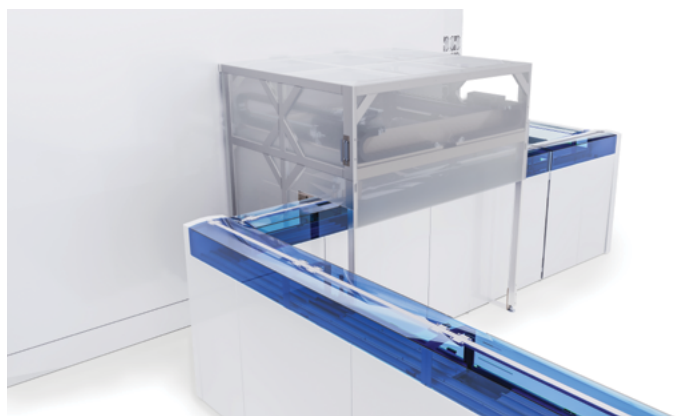

Optional direct connection to TLA is also available

High Throughput

- Meets the needs of the largest labs in the world
- Loads up to 10,000 tubes per hour
- Automated high-throughput discard

Powerful Software

Verso Dx offers automatic sample tracking and access controls for every processed specimen with our proprietary INSTINCT S automation software. A complete audit trail is maintained, and various reports can be generated with a few clicks. While the software is incredibly powerful, daily tasks are easy to learn and lab personnel can operate the system with minimal training.

Additionally, the RESTful LIS interface offers enhanced connectivity to LIS systems.

Compact Footprint

- Smaller than other automated systems, manual freezer farms, and walk-in fridges
- Flexible system heights available



High capacity for manual input: up to 1,300 samples can be loaded at once



● Headquarters / Manufacturing



Years of Experience
75+



Locations Worldwide
22+



Employees Internationally
3,000+

To find a representative in
your area, please visit:
www.hamiltoncompany.com/contacts

Web: www.hamiltoncompany.com
USA: 800-310-5866
Europe: +41-58-610-10-30

Americas & Pacific Rim
Hamilton Storage
Technologies, Inc.
3 Forge Parkway
Franklin, MA 02038 USA
Tel: +1-508-544-7000
sales@hamilton-storage.com

Europe, Asia & Africa
Hamilton Storage GmbH
Parc Industrial Vial 10
7013 Domat/Ems, Switzerland
Tel: +41-58-610-10-30
sales@hamilton-storage.com

Japan
Hamilton Company Japan
Shiroyama Trust Tower 31F
4-3-1 Toranomon, Minato-ku
Tokyo, 105-6031, Japan
Tel: +81-3-6435-6850
sales@hamilton-storage.com

